SIEMENS

Transmitter

User Guide


usa.bestsoundtechnology.com

Life sounds brilliant.

Enhance your Siemens experience.

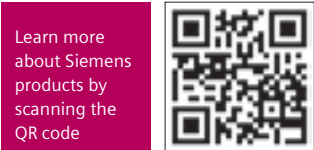Learn about your hearing aids and accessories you may find helpful, access the latest information about hearing health, and receive promotional offers by joining the Siemens hearing health community at* usa.siemens.com/join or:


*We respect your privacy. We will not sell, rent, or give your name or e-mail address to anyone. At any time, you can select a link at the bottom of every email to unsubscribe or update profile.

Content

Transmitter	2
Included in delivery	2
Intended use	3
General notes on safety	3
Basic functions	6
Connecting Transmitter to easyTek™	8
Connecting Transmitter to miniTek®	9
Connection Scheme	11
Connecting the Transmitter to your television	12
Maintenance and care	13
Disposal information	13
Bluetooth®	14
Technical information	14
Technical Data	16

Symbols:	
Be aware of information marked with the warning symbol WARNING , CAUTION or NOTICE !	
	WARNING points out a situation that could lead to serious injuries, CAUTION indicates a situation that could lead to minor and moderate injuries.
	NOTICE indicates possible property damage.
	Advice and tips on how to handle your device better.
	Instruction. Indicates that something has to be done.

Transmitter

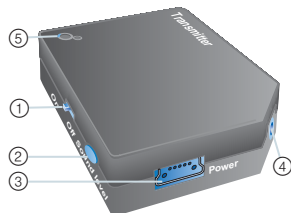- ① On/off switch
- ② Sound level
- ③ USB-socket for power plug
- ④ 2.5 mm stereo socket
- ⑤ LED

Included in delivery

USB-power plug with country specific adapters

Cables:

- 1.8 m audio cable with 3.5 mm stereo plug and a 2.5 mm stereo plug
- 0.3 m audio cable with 3.5 mm stereo plug and a 2.5 mm stereo plug

Adapter set:

- 2.5 mm stereo plug to 3.5 mm stereo socket
- 6.3 mm stereo plug to 3.5 mm stereo socket
- RCA (cinch) to 3.5 mm stereo socket
- RCA (cinch) to SCART

Intended use

Transmitter is intended to make devices without *Bluetooth*® function compatible to *Bluetooth*.

Read and follow the instructions of this user guide to avoid injuries and to protect your Transmitter from damage.

General notes on safety

WARNING

Choking hazard!

Your device contains small parts which can be swallowed.

- Keep hearing instruments, batteries and accessories out of reach of children and mentally disabled persons.
- If parts have been swallowed consult a physician or hospital immediately.

WARNING

Wireless systems may interfere with measuring devices and electronic equipment.

- ▶ Turn your device off when in hospitals or airplanes.
- ▶ Prior to use in conjunction with an electronic implant or life supporting system, have the electromagnetic compatibility verified.

WARNING

Risk of interference with active implants – e.g. pacemakers.

- ▶ Consult a cardiologist before using accessories with audio streaming functions together with active implants.

NOTICE

Your device is sensitive to extreme heat, high humidity, strong magnetic fields (> 0.1T), X-rays and mechanical stress.

- ▶ Do not expose your device to extreme temperature or high humidity.
- ▶ Do not leave it in direct sunlight.
- ▶ Do not immerse it in water.
- ▶ Do not place your device in a microwave oven.

WARNING

Risk of explosion!

- ▶ Do not use your device in explosive atmospheres (e. g. in mining areas).

WARNING

Risk of electrical shock!

- ▶ Do not use obviously damaged devices and return them to point of sale.

WARNING

Risk of explosion when the rechargeable battery is handled wrong.

- ▶ Do not short circuit, pierce, crush or disassemble the rechargeable battery.
- ▶ Prior to return the battery for recycling, tape the electrodes to prevent short circuit, which may cause the battery to explode.
- ▶ Do not place it in oven or fire.
- ▶ Do not immerse it in water.
- ▶ Charge the rechargeable battery between 32 and 104°F.

WARNING

Risk of electric shock!

- ▶ Connect audio input only to a device conforming to IEC 60065, the IEC-Norm for audio, video and other electronic devices.

Basic functions

Before first use

Your Transmitter is equipped with a rechargeable battery.

- ▶ Charge your Transmitter if necessary using the supplied USB-power plug.

LED state while charging:

	Device on	Device off
Audio mode	flashes red-blue	shines red

Charging time is about 3 hours depending on the status of the rechargeable battery.

WARNING

Risk of electric shock!

- ▶ For safety reasons use the supplied USB power supply only.

- ▶ Clean surface before applying clip or adhesive tape.

- ▶ Use the supplied adhesive tape to affix the Transmitter in a suitable place.

Turning on/off

- ▶ Use the on/off switch to turn the Transmitter on or off.

Changing the volume (Sound level)

- ▶ Press the sound level button to change the volume of the audio signal.

The sound level can be adjusted in 5 levels.

The LED flashes blue. The number of times the LED flashes indicates the volume level.

Connecting Transmitter to easyTek™

Using our transmitters

Every transmitter has to be paired once with easyTek. You can use up to two transmitters.

Preconditions for pairing:

- easyTek is turned on. Press  until the green LED turns on . This takes approximately 2 seconds.
- Neckloop is plugged into both sides, easyTek is worn around the neck.

Pairing transmitters:

- ▶ Turn on the transmitter.
- ▶ Press ,  and  simultaneously until the blue and green LEDs flash . This takes approximately 5 seconds.

easyTek searches for transmitters nearby.

Pairing may take several seconds. The blue LED indicates the pairing status:

Pairing status for transmitters

flashing, then on

- ▶ Pairing has been completed successfully. Streaming starts automatically. **You can now use the transmitter.**

flashing, then off

- ▶ Pairing was not successful. In this case, ensure that the transmitter is within operating distance and try again.

You can also use the easyTek App for pairing

Toggle between transmitters:

-  double press for transmitters. Or to toggle transmitters.

Further information:

- If the transmitter is only temporarily (< 5 min) out of operating distance, the transmission of the audio signal is interrupted but automatically picked up again when the transmitter is back within operating distance.
- If a transmitter is already paired, you cannot pair it again.

Connecting Transmitter to miniTek®

- ▶ Connect the Transmitter with the power plug to a power supply or charge it regularly.
- ▶ Use one of the supplied cables to connect the Transmitter to an external device (refer to graphic “Connection scheme”).
- ▶ Press on/off switch to turn the Transmitter on.

The LED on the Transmitter flashes while searching a connection to your streaming device (remote control).

- ▶ Press  for about 3 seconds to turn the remote control on.
- ▶ Press  on your streaming device (easyTek or miniTek) to change to the audio source manually.

 flashes while your streaming device is searching for a connection and it shines constantly if a connection is established.

The LED on the Transmitter shines blue if a connection is established.

Each time you press , your streaming device switches to the next available audio source.


► Start the audio transmission on your streaming device.

You can now hear the sound from your external streaming device through your hearing instruments.


► Press  on your streaming device to stop listening to the external audio source.

Your hearing instruments switch back to the hearing program you used before.

Connection Scheme

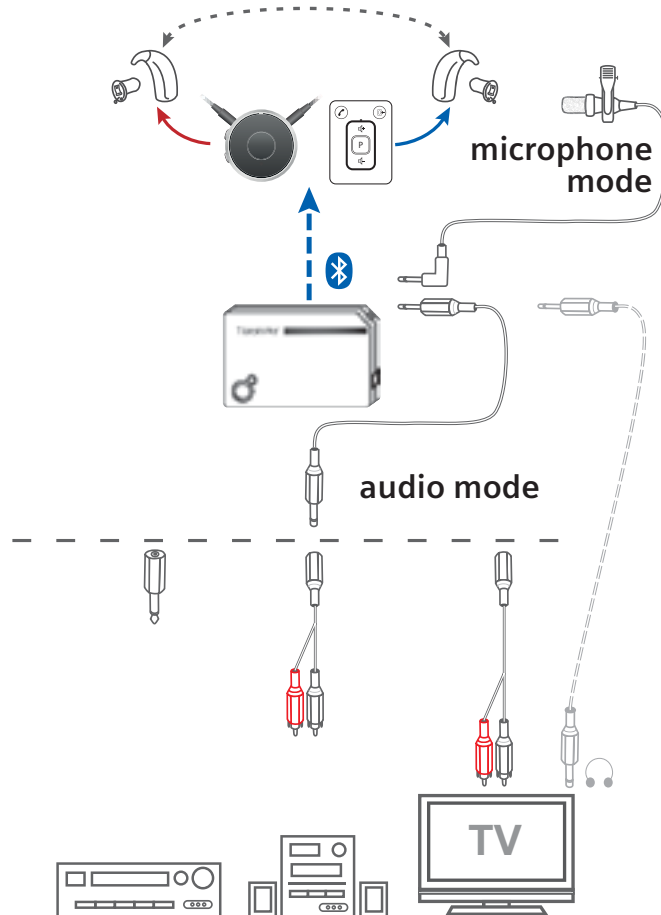Connecting the Transmitter to your television

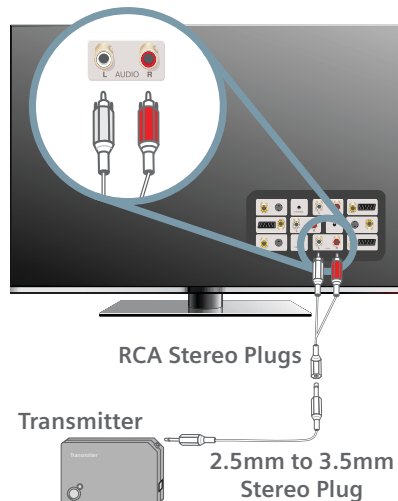Use the RCA adapters with the analog outputs on your TV, cable to satellite dish console, or A/V receiver to connect the transmitter.

Not all television sets support an analog audio output.

- ▶ Ask the dealer or an experienced radio/TV technician for a SPDIF/TOSLINK D/A converter.

If needed, a digital-to-analog converter can be purchased.

Maintenance and care

NOTICE

- ▶ Do not put your device in water!
- ▶ Do not clean your device with alcohol or benzene.
- ▶ Clean the Transmitter as necessary with a soft tissue.
- ▶ Avoid abrasive cleaning liquids with additives such as citric acids, ammonia, etc.

Disposal information

- ▶ Recycle hearing instruments, accessories and packaging according to local regulations.
- ▶ To avoid environmental pollution, do not throw batteries into household trash.
- ▶ Recycle or dispose of batteries according to local regulations or return them to your Hearing Care Professional.

Bluetooth®

Bluetooth is a wireless standard that provides a way to connect wireless devices like cell phones, TV, headsets and others.

The *Bluetooth* connection has a range up to 32 ft (10 m). Walls or objects between the devices can shorten the distance.

► Refer to the user guide of your connected device for more detailed information on how to link devices via *Bluetooth*.

Technical information

Transmitter

FCC: QWOBTA-131

IC: 4460A-BTA131

Federal Communication Commission Interference Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- This device may not cause harmful interference, and
- This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio

frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Caution: Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Industry Canada statement

This device complies with RSS-210 of the Industry Canada Rules. Operation is subject to the following two conditions:

- This device may not cause harmful interference, and
- This device must accept any interference received, including interference that may cause undesired operation.

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Information in this brochure is subject to change without notice. Copyright © 2015 Siemens AG. All rights reserved.
La información que contiene este folleto está sujeta a cambios sin aviso. Copyright © 2015 Siemens AG. Todos los derechos reservados.
under Trademark License of Siemens AG
Sivantos GmbH
Sivantos, Inc.